ARDESTO

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IMPORTANT SAFEGUARDS

Please read the safety instructions before using this product. The safety instructions provided here are to help you use this product safely and correctly, and also to prevent you from injuring yourself or others.

DANGERS-VOLTAGE

- The equipment includes a micro USB cable: length 0.5 m, compatible with any 220 V / 5 V 600 mA adapter. The adapter is not supplied.
- Voltage converters do not guarantee compatibility.
- Do not force the plug into an outlet.
- If the plug does not fit easily into the outlet, discontinue use.
- Always unplug this product immediately after using.

MEDICAL-WARNINGS

- If you have had oral surgery or gum surgery within the past two months, consult your dentist before using the product.
- Consult your dentist if excessive bleeding occurs while using this product or minor bleeding continues to occur after one week of use.
- Consult your physician prior to using the toothbrush if you have any medical concerns.
- Consult your physician or implant device manufacturer prior to using this product if you have questions concerning the combined use of this product with a pacemaker or other implanted medical device(s).

WARNINGS

- This device can be used by children over 8 years of age and persons with physical or sensory or mental disabilities or with insufficient experience and knowledge, if they are under constant supervision or instructed in the safe use of the device and they understand the possible dangers.
- Children should be supervised to ensure they do not play with the product. Cleaning and customer service should not be performed by children without supervision.
- Use this product only for its intended use as described in this manual. This product is designed to clean your teeth, gums only.

Do not use any attachments or brush heads other than those recommended by the manufacturer.

- The toothbrush is a personal care appliance and is not intended for use on multiple patients in a dental office or institution.
- If your toothpaste includes peroxide, baking soda or other bicarbonate (common in whitening toothpastes), thoroughly clean your brush head and handle to avoid cracking the plastic.
- Never operate this product if it has a damaged cord, plug or charger. Do not use this product if it is not working properly, if it has been dropped or damaged or dropped into water. Immediately discontinue use if the product appears damaged in any way.
- Keep the product away from heated surfaces.
- If you experience discomfort or pain in your mouth, discontinue use of this product and contact a dentist.
- Do not attempt to remove the battery. The device has a built-in rechargeable battery that must be disposed of separately only after the end of the life of the electric toothbrush. More details in the section «Disposal».
- Do not clean the product in the dishwasher.
- Never drop or insert any object into any opening.


Production date and device serial number are indicated on the packaging.

SAVE THESE INSTRUCTIONS

EQUIPMENT

- Electric toothbrush - 1 pc.
- Brush head - 1 pc.
- Micro USB cable - 1 pc.
- Operating instructions (contains a warranty card) - 1 pc.

CONFIGURATION

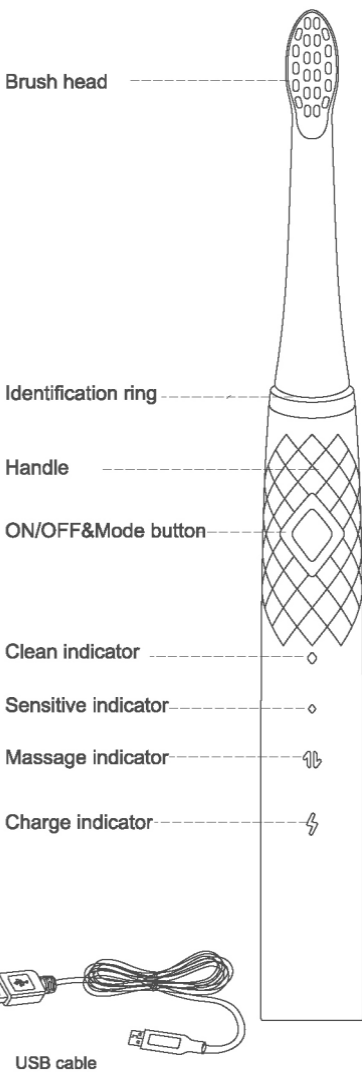GETTING STARTED

Charging the toothbrush

- Be sure to turn off the toothbrush before charging.
- Connect the USB cable to power supply, the charging indicator will flash (Fig. 1).

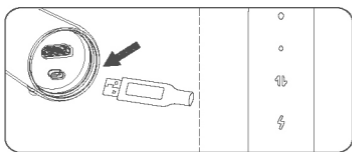
Fig. 1

- It takes at least 3-4 hours to charge the battery fully.
- The toothbrush must be charged for 3-4 hours, if:
 - a) this is first time using the toothbrush;
 - b) the toothbrush has not been used for more than 3 months.
- During charging, the toothbrush will get warm, but this is normal.
- The operating duration will become shorter as the battery gets ageing. If the toothbrush operates only for one short time after being fully charged, this means that the battery is close to the end of its life cycle. Contact service center for support.

HOW TO USE

1. To mount the brush head on the toothbrush, insert it into the shaft. (Fig. 2).

- The shaft prevents the brush head from becoming detached due to vibration.

2. Put toothpaste on the brush head.

3. Place the brush head in your mouth.

4. Light press the button on the unit (Fig. 3), then you can start brushing tooth, at the same time the mode indicator will light.

CLEAN: standard mode for superior teeth cleaning.

SENSITIVE: gently cleans your teeth or any other sensitive areas.

MASSAGE: spulsates for gentle gum stimulation.

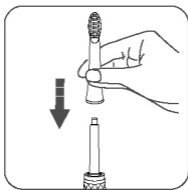
Fig. 2

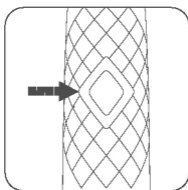
Fig. 3

- Brush evenly throughout the mouth in four sections. (Fig. 4).

1. Outside upper;
2. Outside lower;
3. Inside upper;
4. Inside lower;

5. Guide the brush head slowly from tooth to tooth.

- Hold the brush head in place for a few seconds before moving on to the next tooth.
- Brush the gums as well as the teeth, first the outsides, then the insides, finally the chewing surfaces.
- Do not press too hard or scrub, simply let the brush do all the work.

6. During brushing, the handle will work at the function of 2 minutes cycle with 30 seconds intervals. the unit will switch to next mode, or it will turn off. (Fig. 5 / Fig. 6).

NOTE:

- A brushing time of 2 minutes is generally long enough to clean your teeth properly. If you want you can brush longer than 2 minutes for an even more thorough clean.
- When you first use you electric toothbrush you may experience slight bleeding. This should stop after a few days when your gums get used to the toothbrush. If the bleeding persists, Please consult your dentist.

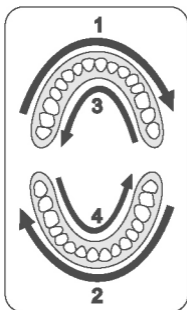
Fig. 4

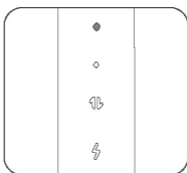
Fig. 5

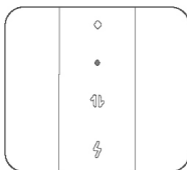
Fig. 6

MAINTENANCE

1. Turn on the toothbrush and rinse the brush head under a tap (Fig. 7).
2. Turn off the toothbrush detach the brush head and rinse the inside of the brush head (max. 140F). (Fig. 8).
3. Turn on the toothbrush and rinse the drive shaft under the tap to remove toothpaste residues. (Fig. 9).
4. Toothbrush head put back on the handle.
5. Clean the toothbrush with a soft moist cloth or paper towel.

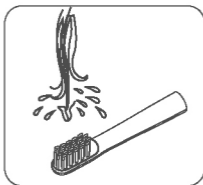
Fig. 7

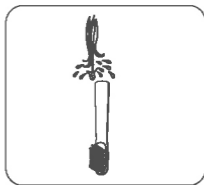
Fig. 8

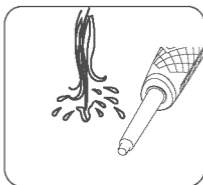
Fig. 9

NOTE:

- Replace the brush head every 3 months for the best brushing effect.
- Also replace the brush head when the bristles start to bend outwards.
- Any other servicing should be performed by an authorized service representative or that the product has no user serviceable parts.

TROUBLESHOOTING

Problem	Possible cause	Solution
The toothbrush doesn't work.	You've just purchased the product or it has been left unused for more than 3 months.	Charge the main unit for at least 3-4 hours.
Mode indicator light, but the toothbrush doesn't work.	Battery is low.	Charge the main unit for at least 3-4 hours.
Weak or diminished oscillations.	The toothbrush is pressed hard against your teeth.	Place the brush head gently against your teeth.
It only works for a few minutes even after charging.	The battery has reached the end of its life cycle.	Contact service center for support.
	Charging time was insufficient.	Charge the main unit for at least 3-4 hours.

DISPOSAL

The power source of this product is Lithium ion hydride battery which can be recycled. At the end of the product's life, always remove the battery before you discard the product, and recycle or dispose of the battery according to local waste management requirements. If you have recycling questions, contact your local waste management office.

Removing and recycling the battery

1. Removing the glue around at the gap between identification ring and handle with a knife (Fig. 10).
2. Bending the shaft in different direction with a little strength unit it comes out of the handle, then cut off the wires with a scissors. (Fig. 11).
3. Clamp and press the handle with a pliers on the position about 3-4 cm from the bottom and push the post of the bracket till it comes out. (Fig. 12)
4. Melting the tin that connects the battery and PCBA with a electric iron, take the battery out and recycle it.

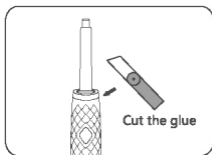
Fig.10

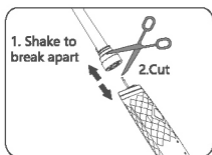
Fig.11

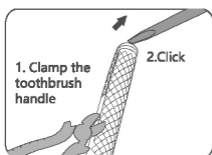
Fig.12

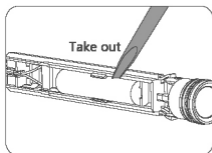
Fig.13

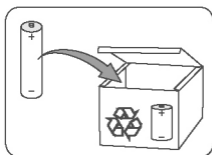
Fig.14

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Warranty conditions:

The plant guarantees normal operation of the product within 12 months from the date of the sale, provided that the consumer complies with the rules of operation and care provided by this Manual. Service life of the product is 3 years.

In order to avoid any misunderstanding we kindly ask you to carefully study the Owner's Manual, the terms of warranty, and check the correctness of filling the Warranty Card. The Warranty Card is valid only if the following information is correctly and clearly stated: model, serial number of the product, date of sale, clear seals of the company-seller, signature of the buyer. The model and serial number of the product must correspond to those specified in the Warranty Card.

If these terms are violated, and if the information specified in the Warranty Card is changed, erased or rewritten, the Warranty Card is considered invalid.

Setup and installation (assembly, connection, etc.) of the appliance are described in the attached documentation; they can be carried out both by the user, and by specialists having the necessary qualification sent by USC or selling companies (on a paid basis). In this case, the person (the institution) installing the appliance shall be responsible for the correctness and quality of installation (setup). Please pay attention to the importance of the correct installation of the appliance, both for its reliable operation, and for obtaining warranty and free service. Require the installation specialist to enter all the necessary information about the installation in the Warranty Card.

In case of failure of the appliance during the warranty period by the fault of the manufacturer, the owner is entitled to free warranty repair upon presentation of the correctly filled Warranty Card together with the appliance in the factory configuration packed in the original packaging to the warranty workshop or the place of its purchase. Satisfaction of the customer's claim due to the fault of the manufacturer shall be made in accordance with the law «On Protection of Consumers' Rights». In the case of warranty repair the warranty period shall be extended for the period of repair and shipment.

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WARRANTY CARD

Dear Buyer! Congratulations on your purchase of the ARDESTO brand appliance, which was designed and manufactured in accordance with the highest quality standards, and we thank you for choosing this particular instrument.

We ask you to keep the coupon during the warranty period. When purchasing a product, require a full warranty card.

1. Warranty service is carried out only if there is a correctly and clearly filled original warranty card, which indicates: product model, date of sale, serial number, warranty service period, and the seller's seal. *
2. The service life of household appliances is 3 years.
3. The product is intended for consumer use. When using the product in commercial activities, the seller / manufacturer does not bear warranty obligations, after-sales service is performed on a paid basis.
4. Warranty repair is carried out within the period specified in the warranty card for the product in an authorized service center on the conditions and terms determined by applicable law.
5. The product is withdrawn from the warranty in case of violation by the consumer of the operating rules set forth in the instruction manual.
6. The product is removed from warranty service in the following cases:
 - misuse and non-consumer use;
 - mechanical damage;
 - damage caused by the ingress of foreign objects, substances, liquids, insects;
 - damage caused by natural disasters (rain, wind, lightning, etc.), fire, domestic factors (excessive humidity, dust, aggressive environment, etc.)
 - damage caused by non-compliance of power and cable network parameters with state standards and other similar factors;
 - when operating equipment in the power supply network with a missing single ground loop;
 - in case of violation of seals installed on the product;
 - lack of serial number of the device, or inability to identify it.
7. The warranty does not cover consumables and accessories, as well as filters, shelves, drawers, stands and other containers for storing products.
8. The warranty period is 12 months from the date of sale.

* Tear-off maintenance tickets are provided by an authorized service center.

The completeness of the product is checked. I have read the terms of the warranty service, no complaints.

Customer Signature _____

Elit service :

3, Petre Kavtaradze str., Tbilisi, 94,

Gorgiladze str., Batumi

71, Davit Guramishvili str., Gori,

73 Kostava str., Akhaltsikhe

Product Information

Product

Model

Serial number

Seller Information

Name of trade organization

Address

Date of sale

Seller stamp

Warranty coupons

Ticket number 3

Saller stamp

date of the application

Cause
of damage

Date
of execution

Ticket number 2

Saller stamp

date of the application

Cause
of damage

Date
of execution

Ticket number 1

Saller stamp

date of the application

Cause
of damage

Date
of execution